



## SpeedScrew® Fully Threaded OPUS® Knotless Fixation Implant

Sophisticated Art for the Shoulder.  
Elite Performance for the Best Outcome.

### Advancing Your Artistic Vision

- The 5.5mm fully threaded design locks securely within the bone to maximize fixation strength
- Helps prevent the possibility of anchor "pull-back" by setting bone lock prior to suture lock
- Suture lock does not depend on bone quality because it is locked within the implant reducing the possibility of suture slip, gap formation or repair failure
- Independent tensioning allows the surgeon to tension suture limbs either individually or together, to evenly distribute pressure on the tendon and to provide maximum lateralization of the tissue while minimizing gap formation
- Radiolucent implant made with PEEK-OPTIMA® polymer is MRI and X-Ray compatible



SpeedScrew  
OM-6500



SpeedScrew Knotless Suture  
Bridge with Doubleplay®

The OPUS System from  
ArthroCare® Sports Medicine



**SpeedScrew**

The latest fully threaded  
implant technology



**Magnum® X**

Independent tensioning for  
the ultimate control



**TwinLock®**

For the most challenging  
rotator cuff repairs



**Magnum PI**


PEEK anchor for Peak  
performance



**Magnum<sup>2</sup>**

The ultimate in speed and  
strength

## SpeedScrew Anchor Specifications

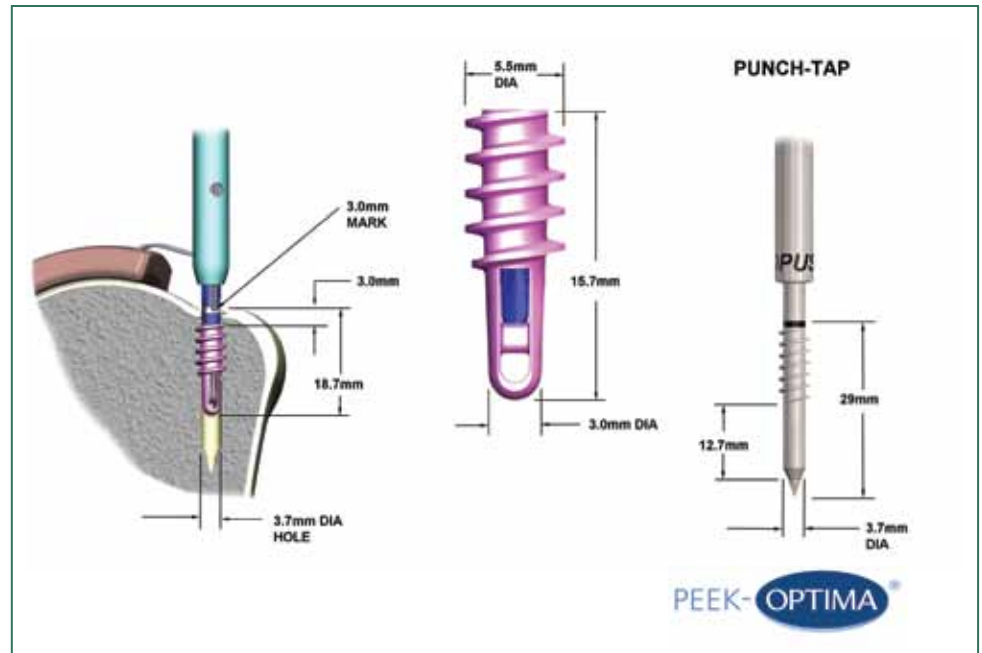

## SpeedScrew Implant Configurations

Part Number	Description
OM-6500CS	5.5mm SpeedScrew Implant, Case
OM-9610	SpeedScrew Insertion Guide
OM-9615	SpeedScrew Obturator/Insertion Guide Set
OM-9620	5.5mm SpeedScrew Punch/Tap
OM-9630	Anchor Extraction Tool
OM-9650	SpeedScrew Obturator
OM-9166	5.5mm SpeedScrew Plus PerfectPasser® MagnumWire® (White) Suture Cartridge
OM-9167	5.5mm SpeedScrew Plus PerfectPasser MagnumWire (Cobraid) Suture Cartridge
OM-9066	5.5mm SpeedScrew Plus M-Connector MagnumWire (White) Suture Cartridge
OM-9067	5.5mm SpeedScrew Plus M-Connector MagnumWire (Cobraid) Suture Cartridge
OM-9161	5.5mm SpeedScrew Procedure Set with PerfectPasser MagnumWire
OM-9061	5.5mm SpeedScrew Procedure Set with M-Connector MagnumWire

To order or for more information, please contact your ArthroCare Sports Medicine representative or customer service at 800-797-6520. Visit our website at [www.arthrocaresportsmedicine.com](http://www.arthrocaresportsmedicine.com).



**ArthroCare Sports Medicine**

680 Vaqueros Avenue  
Sunnyvale, CA 94085-3523  
main phone 408-736-0224  
order entry **800-797-6520**  
order entry fax 888-994-2782  
[arthrocaresportsmedicine.com](http://arthrocaresportsmedicine.com)

**ArthroCare Europe AB**

Skeppsbron 2  
111 30 Stockholm  
Sweden  
phone +46 8 546 172 00  
fax +46 8 546 172 39  
[info@arthrocare.se](mailto:info@arthrocare.se)

**CAUTION:** Federal (USA) law restricts this device to sale by or on the order of a physician.

© 2010 ArthroCare Corporation. All rights reserved. All trademarks, service marks and logos listed herein are registered and unregistered trademarks of ArthroCare Corporation, its affiliates, subsidiaries, or a third party who has licensed its trademarks to ArthroCare.